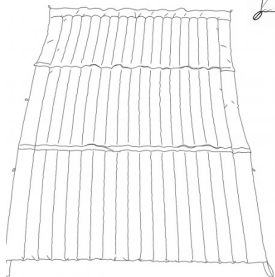
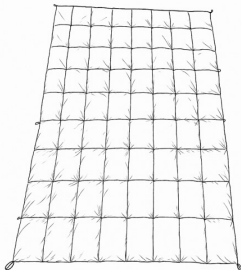
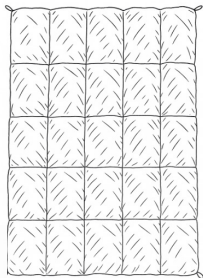
PROTAC®

INTEGRATING SENSES

Protac Kugledynen®

Protac Ball Blanket®

- DA Brugsvejledning
- EN Operations manual
- DE Gebrauchsanweisung
- FI Käyttöohje
- FR Mode d'emploi
- NL Gebruiksaanwijzing
- SV Bruksanvisning
- ES Instrucciones de uso
- IT Istruzioni

| | |
|-----------------|----|
| Dansk | 10 |
| English..... | 15 |
| Deutsch | 20 |
| Suomi | 25 |
| Français | 30 |
| Nederlands..... | 35 |
| Svenska | 40 |
| Español | 45 |
| Italiano..... | 50 |

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English

| | |
|--------------------------------------|----|
| Key to symbols - materials | 15 |
| Marking | 16 |
| Incidents - adverse events | 16 |
| Traceability | 16 |
| Key to symbols - product label | 16 |
| Safety - general precautions | 16 |
| Use of the product | 17 |
| Use environment | 17 |
| Indications | 17 |
| Contra-indications | 17 |
| Maintenance | 18 |
| Cleaning | 18 |
| Key to symbols - care label | 18 |
| Service life and warranty | 18 |
| Disposal | 18 |

You have chosen a Protac quality product. Protac guarantees that this product is of original Danish design and developed in consultation with experts.

Since 1994, Protac has developed, manufactured and sold sensory-stimulating aids to the healthcare sector in Denmark.

Today, Protac products are sold worldwide.

All Protac products have been developed on the basis of the theories of sensory integration, and they are supported by extensive research, documentation and case studies.

The theory of sensory integration by A. Jean Ayres forms the basis of the development of all Protac products. The theory is based on how the brain processes sensory impression and how stimulation of the muscle-joint-sense and the deep tactile sense has a calming and organising effect on the sensory system.

All Protac products contain balls, and it is the movement of the balls that provide the dynamic deep-point pressure on the body and stimulate the receptors in the deep tactile sense and the muscle-joint sense. This enhances the body awareness, alleviates sensory disorder and has a calming effect on both body and mind.

Further information can be found at www.protac.dk.

Key to symbols - materials

100% cotton

100% trevira (flame retardant)

Marking

The products are CE-marked in accordance with Regulation (EU) 2017/745 on medical devices.

EU declaration of conformity can be found at www.protac.dk

Incidents - adverse events

- Incidents related to the use of Protac Ball Blanket® must be reported to Protac.
- Adverse events related to the use of Protac Ball Blanket® must be reported to the National Competent Authority and to Protac.

Traceability

EN Production month, year and batch number can be found on the product label.
Do not remove the product label from the product.

Key to symbols - product label

| | | | |
|---|--|---|-----------------------|
|  | Instructions for use |  | Maximum user weight |
|  | Medical device |  | Article number |
|  | Lot number |  | For indoor use |
|  | Manufacturer |  | Production month/year |
|  | CE marked in accordance with Regulation (EU) 2017/745 on medical devices | | |
|  | Warning: See under Safety | | |

Safety - general precautions

| | |
|---|---|
|  | Blankets with granulate inside contain small particles, and must not be used for children below the age of 3 without supervision. |
|  | The blanket is not to be used as a means of restraint. |
|  | Avoid smoking next to the product. Dropped embers or cigarettes can cause burns. |
|  | Never leave switched on electrical appliances that emit heat on or near the product. |
|  | Do not use near open fire. |
|  | Do not use other balls in the product than those supplied by the manufacturer. |

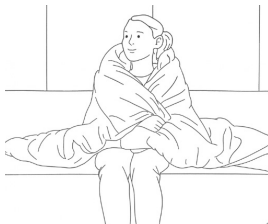
Use of the product

Protac Ball Blanket® is supplied ready for use. It is to be used with a regular duvet cover.

The blanket is intended to be used as a regular blanket for sleeping. The user can lie underneath it or on top of it. For extra sensory stimulation, a helper can "pack" the blanket closer to the body while the user is lying underneath it.

The blanket can also be used for a rest during the day.

EN

Use environment

Protac Ball Blanket® is for indoor use. The blanket can be used both in institutions and private homes.

Indications

Protac Ball Blanket® is a sensory-stimulating aid intended to compensate for a permanent disability.

Protac Ball Blanket® contains balls which move and apply dynamic pressure on to the body, thereby stimulating the muscle-joint sense and the tactile sense. This increases body awareness and calms the brain and body, and the product thereby alleviates physical and mental unrest as well as sleeping difficulties.

Contra-indications

Persons with pressure ulcers or persons at risk of developing pressure ulcers are not to use the product.

Persons unable to care for themselves or persons with a weak lung capacity are not to use the blanket unsupervised.

Maintenance

It is recommended to air the blanket on a regular basis.

Cleaning

Protac Ball Blanket® can be washed in a washing machine. The frequency of washing depends on use. When washing, please follow the instructions on the care label. Washing instructions can also be found on the next page of this manual.

Protac does not recommend the use of stain remover.

EN

Key to symbols - care label

| | | | |
|---|-----------------------------|---|-------------------|
|  | Maximum washing temperature |  | Do not wash |
|  | Can be tumble dried |  | Do not tumble dry |
|  | Dry flat |  | Do not iron |
|  | Line dry in the shade |  | Do not dryclean |
| | |  | Do not bleach |

Service life and warranty

The expected service life of the product is 5 years. This is based on correct use, cleaning, washing and maintenance.

Products used beyond the expected service life must regularly undergo a maintenance inspection which implies checking for wear and tear in fabric and seams.

The product comes with a 2-year warranty, see page 2.

The warranty will be invalid in case of:

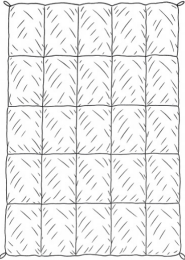
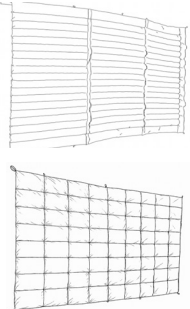
- overload or excessive use of the product
- lack of or incorrect maintenance of the product

Disposal

The product can be broken down into the following components for disposal:

- The plastic balls can be disposed of as hard plastic.
- The fabric can be disposed of as textiles.

| | | | |
|---|--|---|--|
| <p>Protac Ball Blanket® Classic</p> |  | <p>Wash separately at maximum 60°C in a washing machine with a minimum capacity of 16 kg. Use detergent without perfume and no fabric softener. Can be tumble dried at maximum 60°C.</p> <p>Do not dryclean. Do not iron. Do not bleach.</p> | |
| <p>Protac Ball Blanket® Flexible Protac Ball Blanket® Combi</p> |  | <p>Ball bags: Wash separately at maximum 60°C in a washing machine with a minimum capacity of 6 kg. The included laundry bag must be used. Use detergent without perfume and no fabric softener. Tumble dry at maximum 60°C or air dry. Do not dry clean. Do not iron. Do not bleach.</p> | <p>Cover: Wash separately at maximum 60°C in a washing machine with a minimum capacity of 7 kg. Use detergent without perfume and no fabric softener. Can be tumble dried at max. 60°C. Do not dryclean. Do not iron. Do not bleach.</p> |
| <p>Protac Ball Blanket® Calm Protac Ball Blanket® Granulate</p> |  | <p>Unzip the blanket into 3 parts. Wash each part separately at maximum 60°C in a washing machine with a minimum capacity of 7 kg. Use detergent without perfume and no fabric softener. The included laundry bag must be used. Spin at low speed. Can be tumble dried at max. 60°C for max. 10 min., then air dry.</p> <p>Do not dryclean. Do not iron. Do not bleach.</p> | |

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Protac A/S